

Bone Tack

Bone Tack


This Bone tack is designed to be driven in using a mallet at the time of insertion while being screwed out using a driver at the time of removal, with a view to fixing membranes or foils.



BTKIT

- This Bone Tack has a screw structure directly under its head, which allows for easy removal using the driver supplied with the product.
- It is made of titanium GR5 ELI and is available in 3mm and 5mm in size. 3mm Tacks are for normal bone while 5mm Tacks for soft bone.

Consist of Bone Tack

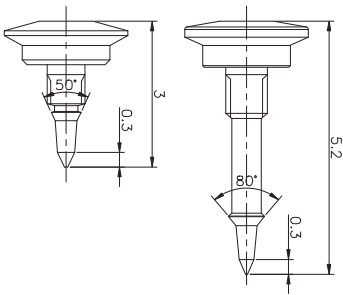

- 1 Bone Tack Holder
- 2 Case
- 3 Bone Tack M0,85 x 5mm
Bone Tack M0,75 x 3mm
- 4 Hexa Driver 0,9mm

Consist of product

		(ea)
Bone Tack Holder	BTSHC(Straight)	1
Case	BTSC	1
Bone Tack M0,85×5mm	BTS85-50	5
Bone Tack M0,75×3mm	BTS75-30	9
Hexa Driver 0,9mm	BTIDH09	1

Bone Tack Offset Holder