

## Patient Return Electrodes

Single Use



Latex Free

Latex Free. All Valleylab products are latex free.

DB 15 ft cord



E7506

50 units/case

### Standard Adult PolyHesive™ Patient Return Electrode

- 2.7 m (9 ft) cord



E7507

50 units/case

### Adult REM PolyHesive™ II Patient Return Electrode

- For patient > 30 lbs (> 13.6 kg)
- 2.7 m (9 ft) cord

## Patient Return Electrodes

Single Use



E7507DB

50 units/case

### Adult REM PolyHesive™ II Patient Return Electrode

- For patient > 30 lbs (> 13.6 kg)
- 4.6 m (15 ft) cord



E7508

100 units/case  
Individually Packaged

### Adult REM PolyHesive™ II Patient Return Electrode

- For patient > 30 lbs (> 13.6 kg)
- Requires E0560 reusable cord

## Patient Return Electrodes

Single Use



E7509

50 units/case  
5 per package  
10 ReSealable packages per case

### Adult REM PolyHesive™ II Patient Return Electrode

- For patient > 30 lbs (> 13.6 kg)
- Requires E0560 reusable cord



E7509B

600 units/case  
50 per package  
12 packages per case

### Adult REM PolyHesive™ II Patient Return Electrode

- For patient > 30 lbs (> 13.6 kg)
- Requires E0560 reusable cord

## Patient Return Electrodes

Single Use



E7510-25

25 units/case

### Infant REM PolyHesive™ II Patient Return Electrode

- For patients 6 - 30 lbs (2.72 - 13.6 kg)
- 2.7 m (9 ft) cord



E7510-25DB

25 units/case

### Infant REM PolyHesive™ II Patient Return Electrode

- For patients 6 - 30 lbs (2.72 - 13.6 kg)
- 4.6 m (15 ft) cord




E7512








12 units/case

### Neonatal REM PolyHesive™ II Patient Return Electrode

- For patients < 6 lbs (< 2.72 kg)
- 2.7 m (9 ft) cord

10045586	Puritan Bennett™ cylinder mount bracket for Puritan Bennett™ 840 compressor cart	
10001114	Puritan Bennett™ wall air water trap kit for Puritan Bennett™ 840 pole/compressor cart	
10081875	Hudson™* and Neptune™* humidifier mount kit	
4-076405-00	Neonatal filter adapter (Reusable)	
4-076408-00	DAR™ neonatal D-X800 expiratory filter disposable (12 pack). Single use only.	
4-070305-00	Puritan Bennett™ expiratory bacterial filter (requires drainage vial) (Reusable)	

## PURITAN BENNETT™ 980 and 840 VENTILATOR COMMON PARTS

Part Number	Description	Image
4-MR850-00	Humidifier base; MR850	
4-074600-00	Puritan Bennett™ inspiratory bacterial filter, ISO 22 mm connectors. Single use only.	
4-074601-00	Puritan Bennett™ inspiratory bacterial filter, ISO 22 mm connectors (12 pack). Single use only.	
355U5427	DART™ pediatric-neonatal electrostatic filter HME, small (V <sub>T</sub> 30-100 mL). Single use only.	
10087137	Drain bag holder metal peg	<i>Image not available</i>
4-001474-00	Oxygen hose assembly U.S.	
4-006541-00	Air hose assembly U.S.	
4-048491-00	Drain bag, disposable pack 25. Single use only.	

### 8.15.8.17.2 Installing oxygen and air flow sensors

- 1 Make sure O-rings are properly seated at each end of flow sensor tube (Figure 8-59).
- 2 Orienting flow sensor body as shown in Figure 8-58, gently push sensor into PSOL manifold.
- 3 Install PSOL manifold (Section 8.15.8.16.2).
- 4 Applying leak detector fluid (P/N 4-004489-00) to bottom of flow sensor only, verify that sensor does not leak.

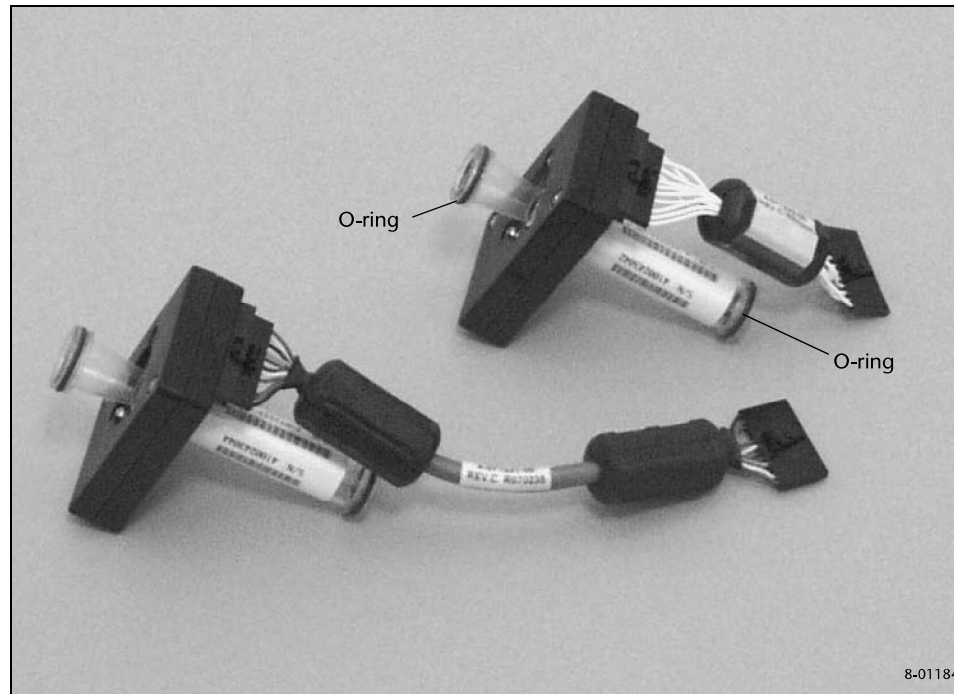


Figure 8-59. Oxygen and air flow sensors (Q1 and Q2)

### 8.15.8.18 Safety valve

#### 8.15.8.18.1 Removing safety valve (SV)

- 1 Remove oxygen and air flow sensors (Section 8.15.8.17.1).
- 2 Remove two flat-head screws that secure safety valve (Figure 8-58) to PSOL manifold. Gently pull safety valve (Figure 8-60) out of manifold.