


SCISSORS

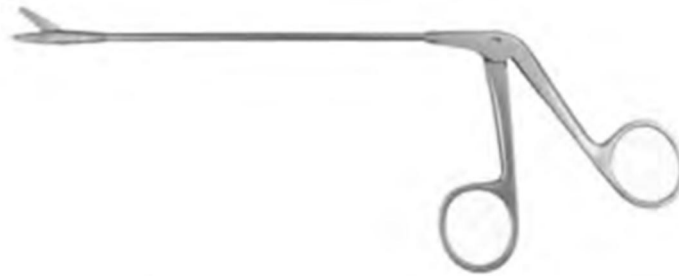

LC-007-001

Other Names: Surgical scissors, suture scissors, nurse's scissors.

Category: Cutting and Dissecting Uses: Used for general cutting of sutures, gauze and other materials; sharp blades are used to cut tougher tissue and are frequently used for dissection; blunt blades are used for cutting materials such as tubing or gauze.

Description: Multi-purpose cutting and dissecting scissors with sharp/blunt, sharp/sharp, or blunt/blunt tip configurations; straight or curved and varying sharpness and lengths.

Notes: The screw must be securely tightened; check for signs of wear and sharpen blades regularly as needed.