



# PeriAcryl®90

HIGH VISCOSITY

Cyanoacrylate Oral Adhesive Rx ONLY

5 mL Multi-Use Kit  
with 50 disposable pipettes and dispensing tray



## PeriAcryl®90 and PeriAcryl®90 High Viscosity Cyanoacrylate dressings that offer surgical predictability and efficiency.

### What are PeriAcryl®90 and PeriAcryl®90-HV?

These medical devices are a special blend of butyl and octyl cyanoacrylate adhesives, formulated to perform in the oral cavity. PeriAcryl®90-HV is the latest dental innovation from GluStitch Inc and has a slower rate of flow. This offers the clinician more control over application. Both products are manufactured in Canada, and North American periodontists have been using PeriAcryl®90 products since 1998.

### Features and benefits of PeriAcryl®90 and PeriAcryl®90-HV:

Coloured in deep violet to increase visibility during application.

Both products polymerize rapidly when they come in contact with water or saliva.

Stays in place for up to 14 days when a multi-layer application is made.

Available in a multi - use kit, specifically designed for dental specialists, to reduce your cost per use. Each kit contains 5 mL of adhesive and 50 disposable pipettes.

PeriAcryl®90 High Viscosity is 9 times thicker than the original formulation which allows for more control during application.

### What are PeriAcryl®90 and PeriAcryl®90-HV used for?

North American periodontists use PeriAcryl®90 products in a variety of surgical situations. These products are effective when used as a liquid bandage to cover free gingival graft donor sites; to help secure grafted tissue; to reinforce sutures; protect sutured sites; and to attach dental dressings.

### Why use PeriAcryl®90 products?

- Simplify surgical procedures
- Increase procedural efficiency
- Improve aesthetic outcomes
- Protect against patient noncompliance
- Increase patient comfort

GluStitch Inc is proudly Canadian owned and operated.

[WWW.GLUSTITCH.COM/PERIACRYL](http://WWW.GLUSTITCH.COM/PERIACRYL)





Search "PeriAcryl" or "GluStitch" in Facebook and you will find our page where we are posting pictures and videos of how periodontists are working with PeriAcryl and the great results they are seeing



email: [info@glustitch.com](mailto:info@glustitch.com) • phone: 1-877-940-2262 • fax: 1-877-450-4000

# How do I use PeriAcryl®90 and PeriAcryl®90-HV?

Both products work the same way. After removing PeriAcryl®90 from the bottle with one of the pipettes supplied, apply a conservative amount overtop of the surgical site and touch it briefly with a wet cotton gauze. In seconds, it will set to form a strong protective covering. In the images below, PeriAcryl®90-HV is the final step in this socket preservation protocol, applied over the membrane and sutures.

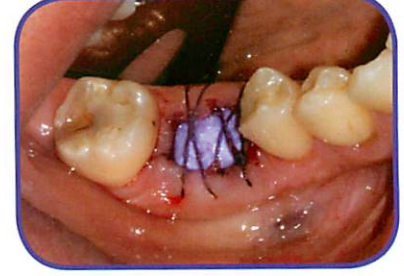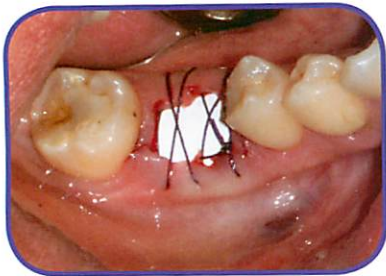

\*donor site 1 week post-op




\*different patient, donor site, 1 week post-op



\*lower right, day of surgery

\*same patient, 1 week post-op

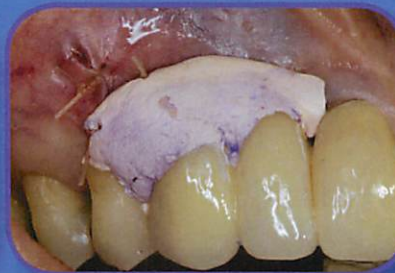

## Where else can PeriAcryl®90 and PeriAcryl®90-HV be used?

Over the donor site of a free gingival graft once hemostasis has been achieved.

Over melanotic lesions post surgery

Over free gingival grafts to reinforce sutures and help secure graft tissue.

Use PeriAcryl®90 or PeriAcryl®90-HV to attach dental dressings like BloodSTOP® and COE-PAK™.