



PARAMOUNT BED


KQ-781

Patient Lifter

**Paramount Bed strives for patient's safety and easy usage.
Our Lifter is cordless for easy operation.**



Hanger hook
Safety stopper is equipped.

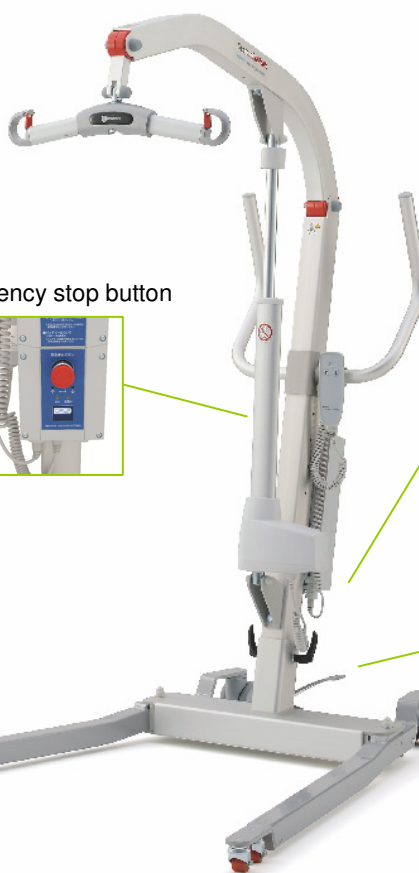


Handle
Height adjustable
5 levels in
increments of 6cm

Height of column is
adjustable at 3 levels
in increments of 5cm.

Emergency stop button



Pedal
Easy to open/close the
legs

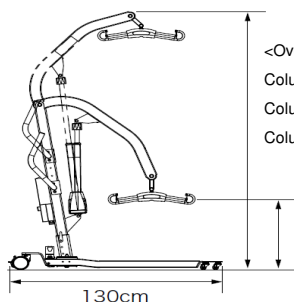

Height of leg is 6.6cm to accommodate to low height bed. Double casters for better maneuverability.

- Material: Column/boom: Aluminum, Base frame: Steel with powder coating
- Safe working load: 205kg
- Weight of product: Approx.38kg
- Dimension: length:130cm / width: 65-106cm



Raise back section, and put sling sheet on patient properly.



Lift up patient by using hand held controller.



Lower patient by using hand held controller, and take off the sling sheet.

*Sling sheet is an option. Please select a sling sheet that suits patient's condition (see back page).



PARAMOUNT BED

KQ-77,78,79A,79B

Sling sheet

Sling seat eases Lift transfer by supporting patient's posture.

KQ-77




Half size

Suitable patient's height: 148-175cm

For patient whose neck is stable.

- polyester mesh type
- Inclusive of a storing bag which can be also used as a washing net.

KQ-78

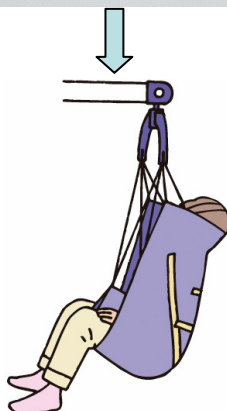


Full size

Suitable patient's height: 148-175cm

For patient whose neck is unstable.

- polyester mesh type
- Inclusive of a storing bag which can be also used as a washing net.

KQ-79A,B

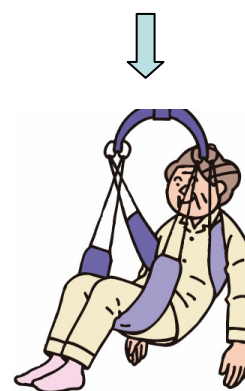

Separate type

Suitable patient's chest:
KQ-79A 60-85cm
KQ-79B 85-100cm

For patient whose head and shoulder are stable.

- Chest/leg sling are sold by one set.
- It is not available for patient whose posture can't be kept when his/her back and leg are lifted.
- polyester mesh type
- Inclusive of a storing bag which can be also used as a washing net.



Scales

Charder MHS2500 is used together with stationary or mobile lifts for weighing of users.

Art. No.	Product name	Description	Unit
70100002	Charder MHS2500	Max 300 kg/ 660 lbs	1 Pce
70100003	Charder MHS2500	Max 400 kg/ 880 lbs	1 Pce

